

Scope

This document describes how to use Health Link software to test PHDC device.

Preparation

Host

Personal computer running Windows Xp or Windows 7.


Device

A board, i.e. twrk70f120m, which is running dev_phdc_weightscale example.

Steps

Follow the steps to install and use Health Link software.

1. Download Health Link software at <http://www.lnihealth.com>
Contact with support@lnihealth.com to get license.
2. Install Health Link software and USB device driver.
3. At the first time run Health Link, you should create new account to store data.
In **Edit** Menu, choose **Advance** -> **Transport Settings** and make sure that **Enable USB Connection** check box is checked

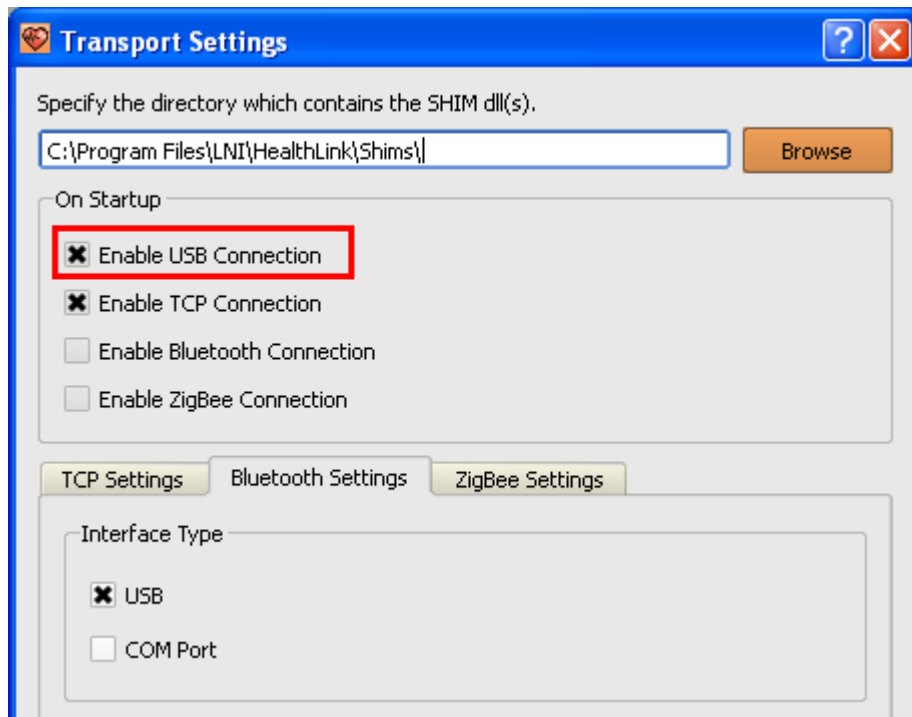

4. Plug-in PHDC device into PC, you can see state of device and data received

